

Human Tau-410 / 2N3R Pre-formed Fibrils Protein, His Tag

Catalog # TAU-H5149



Synonyms

DDPAC, FTDP-17, MAPT, MSTD, MTBT1, Tau, PHF-tau, TAU

Source

Human Tau-410 / 2N3R Pre-formed Fibrils Protein, His Tag (TAU-H5149) is expressed from *E. coli* cells. It contains AA Met 1 - Leu 410 (Accession # [P10636-5](#)).

Predicted N-terminus: Met

Molecular Characterization

Poly-his	Tau(Met 1 - Leu 410) P10636-5
----------	----------------------------------

This protein carries a polyhistidine tag at the N-terminus.

The protein has a calculated MW of 44.5 kDa.

Application

1. Sonication Conditions, Dissolution, Aliquoting, Storage, and Notes for PFF

Thawing: Thaw PFFs rapidly in a 37 °C water bath, or allow to thaw at room temperature.

Aliquoting: Since PFFs are supplied as a suspension, pipette up and down thoroughly before aliquoting to ensure homogeneity.

Storage: Store at -80 °C at all times; avoid storage at 4 °C or -20 °C, which can induce fibril depolymerization. * α -Syn fibrils cold-denatured to monomers at 0-20 °C and heat-denatured at 60-110 °C.

Sonication: The sonication protocols below are based on cell-based assay conditions.

-Option 1 (Probe Sonicator): Use a probe sonicator (SCIENZ) at 10% power (-95 W), applying 60 pulses of 0.5 s on/0.5 s off (recommended).

-Option 2 (Ultrasonic Bath): Sonicate in an ultrasonic cleaner (40 kHz, 200-400 W) at 37 °C for 1 hour; avoid performing the treatment at 4 °C or 20 °C.

Note: The above sonication protocols are based on cell-based assays; other applications (e.g., in vivo injections) may require optimization-users should determine their optimal settings for their specific use case.

2. Recommendations for PFF Use in Animal Models

Minimize Freeze-Thaw Cycles: For animal studies, it is recommend to use the PFF at once or aliquot before use to avoid repeated freeze-thaw cycles.

PFF Maintenance During Injections: During stereotaxic injections, keep the sonicated PFFs in a 37 °C water bath to prevent re-aggregation or sedimentation.

Purity

>90% as determined by SDS-PAGE.

Formulation

Supplied as 0.2 μ m filtered solution in PBS, pH7.4 with trehalose as protectant.

Contact us for customized product form or formulation.

Shipping and Storage

This product is supplied and shipped on dry ice.

Please avoid repeated freeze-thaw cycles.


This product is stable after storage at:

- The product should be stored at -70°C or room temperature for short storage. Do not store fibrils on ice or at 4°C;
- The unsonicated fibril is validated to be stable after storage at -70°C for 1 year under sterile conditions;
- The sonicated fibril should be stored at -70°C for not more than 8 weeks.

ACRO Quality Management System


- [QMS\(ISO, GMP\)](#).
- [Quality Advantages](#)
- [Quality Control Process](#)

SDS-PAGE

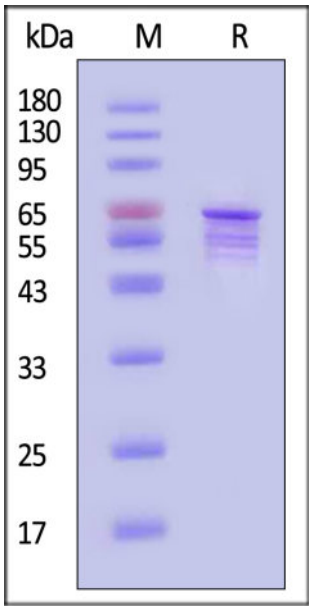

Tau-410 / 2N3R monomer on SDS-PAGE under reducing (R) condition. The gel was stained with Coomassie Blue. The purity of the protein is greater than 90% (With [Star Ribbon Pre-stained Protein Marker](#)).

Electron Microscope

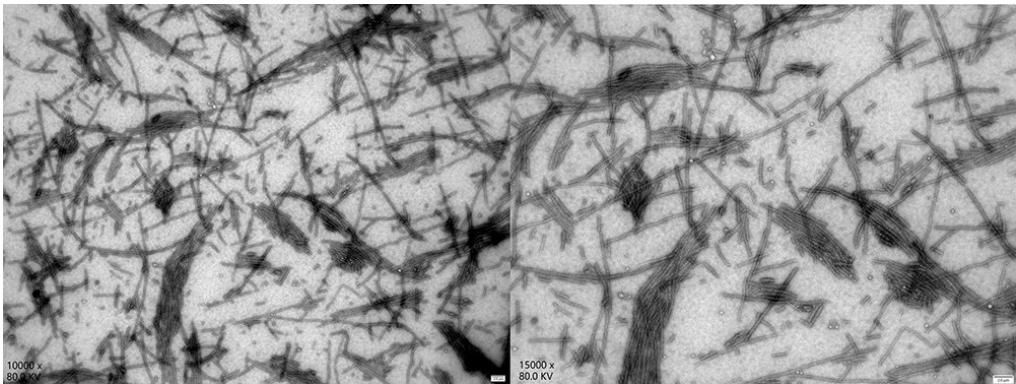

Transmission electron microscopy (TEM) of Tau-410 preformed fibrils (Cat. No. TAU-H5149). Fibrillar structures are visible in negative-stain TEM images of Tau-410 PFFs.

Background

Tau is a microtubule-associated protein, which encodes by the MAPT gene that located on chromosome 17q21. Tau Promotes microtubule assembly and stability, and might be involved in the establishment and maintenance of neuronal polarity. Hyperphosphorylation of the tau protein (tau inclusions, pTau) can result in the self-assembly of tangles of paired helical filaments and straight filaments, which are involved in the pathogenesis of Alzheimer's disease, frontotemporal dementia, and other tauopathies.

Discounts, Gifts,
and more!



Acro



www.acrobiosystems.com